LIBEREX

User Manual


Portable Oral Irrigator

Model NO.: FC 2660

INSTRUCTION (EN) >>

ANLEITUNG (DE) >>


INSTRUCTION (FR) >>

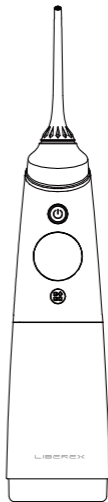Contents

Contents	01
Product Description	02
Product Accessories	03
Operating Procedures	03
Operating Step	05
Recommended Techniques	06
Disassembly Example	07
Battery Maintenance And Handling	07
UL Standard	08
Comment Problems and solutions	10
After-sales Services	11

Product Description

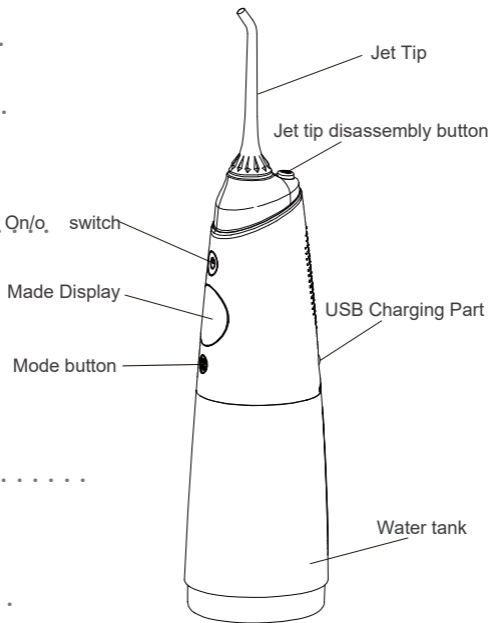Product Accessories

Nozzle x5


Charging Cable x1

Replacement Nozzle / Accessories can be purchased from
LIBEREX on Amazon

Operating Procedures


Usage instruction

Please ensure the machine is full of battery before first time usage, do not use it while charging.

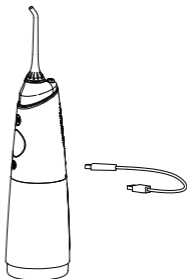Installation and Disassembly of Jet Tip

Press the release button and insert the nozzle into spray head hold until you hear a "click" ring.

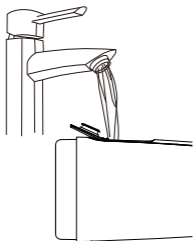Open the water tank lid then fill with the water in it.
Note: Make sure that the device is switched off.

After using, press the release button and put the nozzle.

Operating Step

- ◆ Press the mode button continuously, the mode switches among the sequence: high-medium-low-pulse; This device has a memory function. So in the off state, you can also select the mode suitable for the user and then turn it on again.
- ◆ Irrespective of the mode, after operating for 120 seconds the product will automatically switch off.
- ◆ When power is insufficient, the indicated light will flash until power off. During charging, this product cannot be used.
- ◆ The charger is suitable for: AC 100-240VAC, 50/60Hz, output 5V --- 1A
- ◆ Keep the device vertical for the best effect, rotate the jet tip dial in order to adjust the direction of the water flow.
- ◆ It's normal phenomenon that caused gum bleeding for the first time use. It is recommended to start from a lower pressure level.


Recommended Techniques

Holding the device at downward, aim the jet tip to the teeth by rotating the tip base before switching on the device.


When the water flow and the gums are at a 90 degree angle, slightly close your mouth in order to avoid the water spraying out. In order to achieve the best results. Please start from molar area to the incisors. Slide the water flow along the gums, as well as stay a while between teeth, until the area around the gums and the spaces between the teeth are fully cleaned.

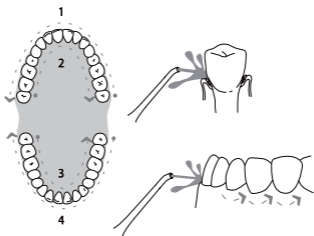
Disassembly Example

Following Usage

Please switch off the power, as well as empty the water tank.

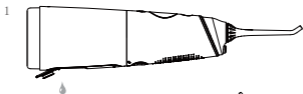Using hands as indicated in Diagram 2 use your left hand to hold the product body, and your right hand to hold the water tank. At the same time rotate downwards in an anti-clockwise direction, in order to remove the water tank.

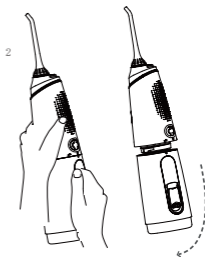Battery Maintenance And Handling

Battery Maintenance

Prior to initial usage the device must to be fully charged, for approximately 4 hours.If to be store for a long-period of time, please fully charge the device in advance.In order to maximize the usage life of the battery, please promptly recharge when power is insufficient.

Battery Handling

The device includes a lithium battery which cannot be changed. At the end of the battery life, please dispose in a local recycling centre.

UL standard

IMPORTANT SAFEGUARDS:

When using electrical products, especially when children are present, the basic safety precautions below should always be followed.

READ ALL INSTRUCTIONS BEFORE USING.

DANGERS:

To reduce the risk of electrocution:

- ◆ Do not handle charger with wet hands.
- ◆ Do not immerse in water or other liquid.
- ◆ Check the charger cord for damage before use.
- ◆ Keep charger and cord away from heated surfaces.
- ◆ Use this product only as indicated in these instructions or as recommended by your dental professional.
- ◆ Do not reach for a product that has fallen into water. Unplug immediately.
- ◆ Do not plug this device into a voltage system which is different from specific voltage range.

To reduce the risk of burns, electric shock, fire, serious injury or damage issues:

- ◆ Do not use if you have trauma on your tongue or in your mouth.
- ◆ If you are in dental treatment, you should consult your dentist before using this instrument.
- ◆ Do not direct water flow under the tongue, into the ear, nose or other delicate organs.
- ◆ It's normal phenomenon that caused gum bleeding for the first time use. It is recommended to start from a lower pressure level.
- ◆ Never be left unattended when plugged in. Close supervision is necessary when it used by, on, or near children or invalids.

Comment Problems and solutions

Malfunction	Cause	Handling method
Water leakage between the jet tip and screw part.	1.The jet tip cannot be installed properly. 2.The seal ring socket for the jet tip is damaged.	1.Insert the jet tip into the screw part until you hear a clicking sound. 2.Contact the retailer for repair and install the jet tip properly
Inadequate pressure	Seal ring is loose or damaged	Replace the seal ring
Inadequate pressure	Low voltage usage	Use after charging
Device cannot be switch on	Low voltage usage	Use after charging
Device is not spray water	Filtration net is clogged	Wet and clean the Filtration net
Water tank sides are leaking water	Air entry hole is producing water	Standard phenomenon high air pressure principle

If the above methods cannot resolve these problems, please consult with the retailer or proceed to the designated after-sales repair site for repairs.

Warranty

Register your product with us in order to get extra 12 months warranty (24 months in total): including free product replacement, technique support and any important updates or upgrades to your product:

1. Please visit our website: frigif.com/liberex
2. Fill in the online form to register and get extra 12 months warranty

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